



Trīs komponentu PHB maisījuma hromatogramma

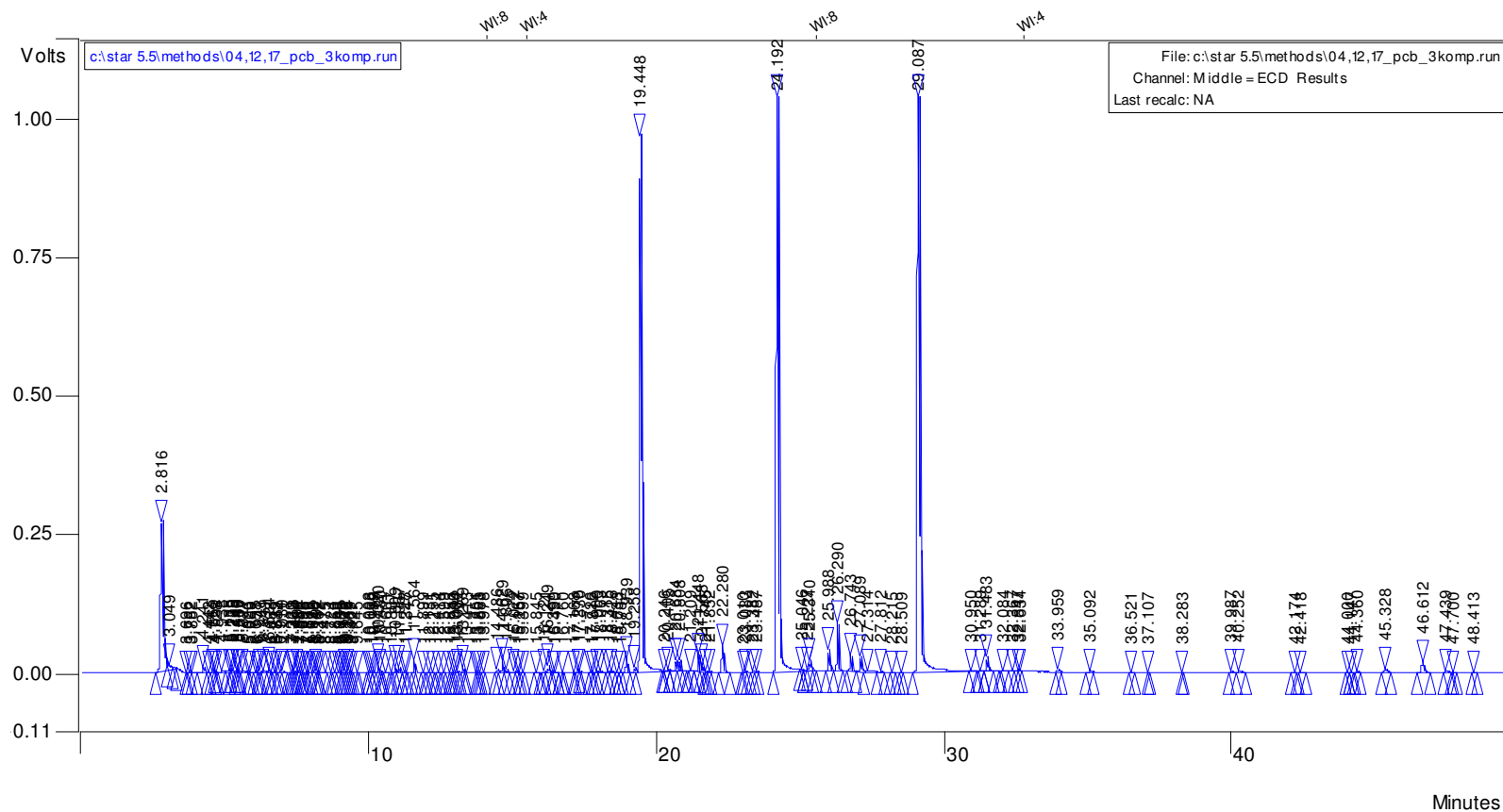

3 PHB Congener Mix C-SCA-01, AccuStandard Vielu izdalīšanas laiki: PHB 77 – 19.46 min, PHB 126 – 24.19 min un PHB 169 – 29.09 min.


Sešu komponentu PHB maisījuma hromatogramma

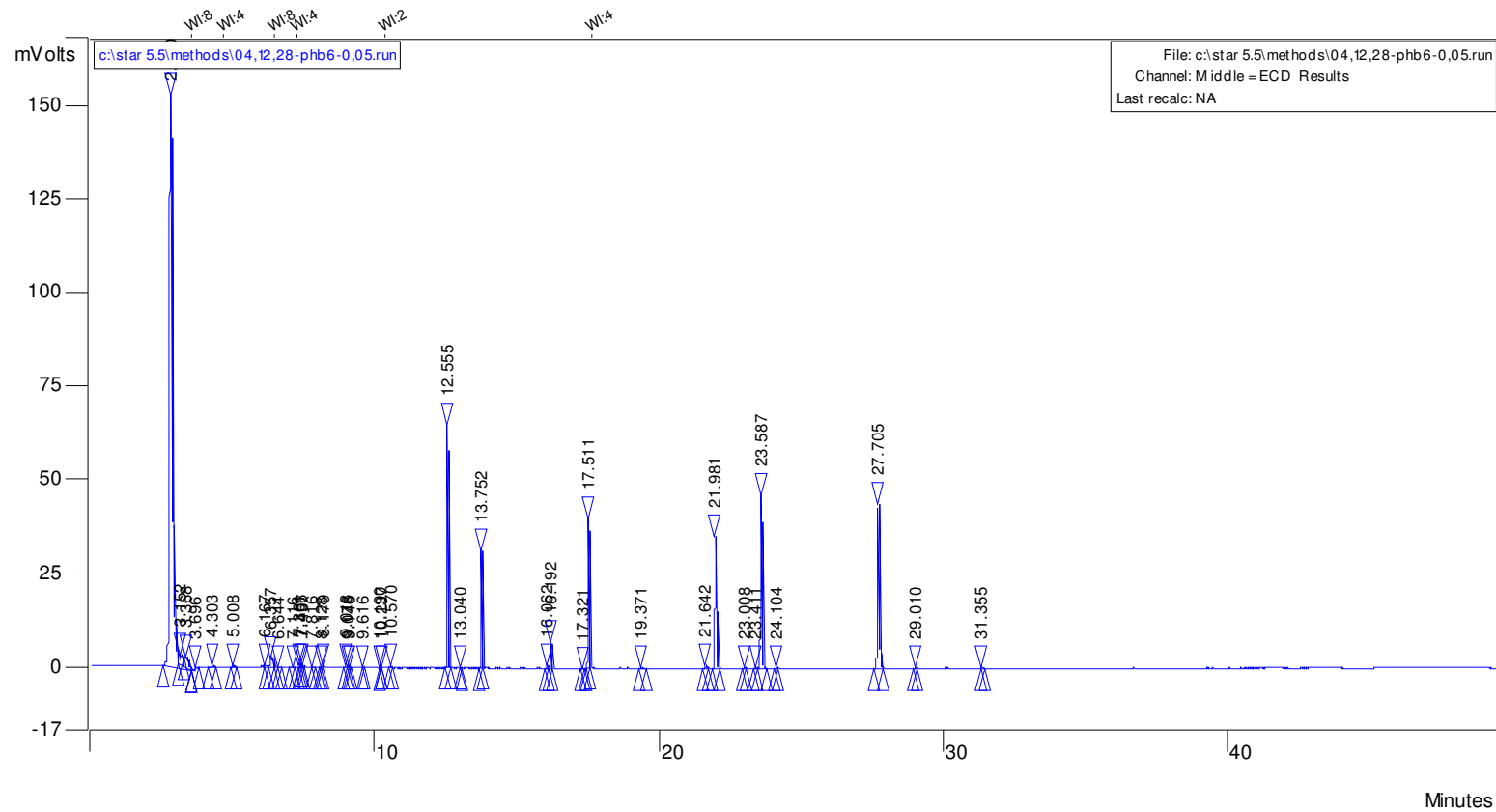

Sešu PHB standarts – PCB Congener Content Evaluation Mix AE-00059, AccuStandard. Vielu izdalīšanas laiki: PHB28 – 12.56 min, PHB 52 -13.8 min, PHB 101 – 17.51 min, PHB1153 – 22.0min, PHB138 – 23.6 min, PHB 180– 27.71 min.

8. pielikums

Divpadsmit komponentu PHB maisījuma hromatogramma

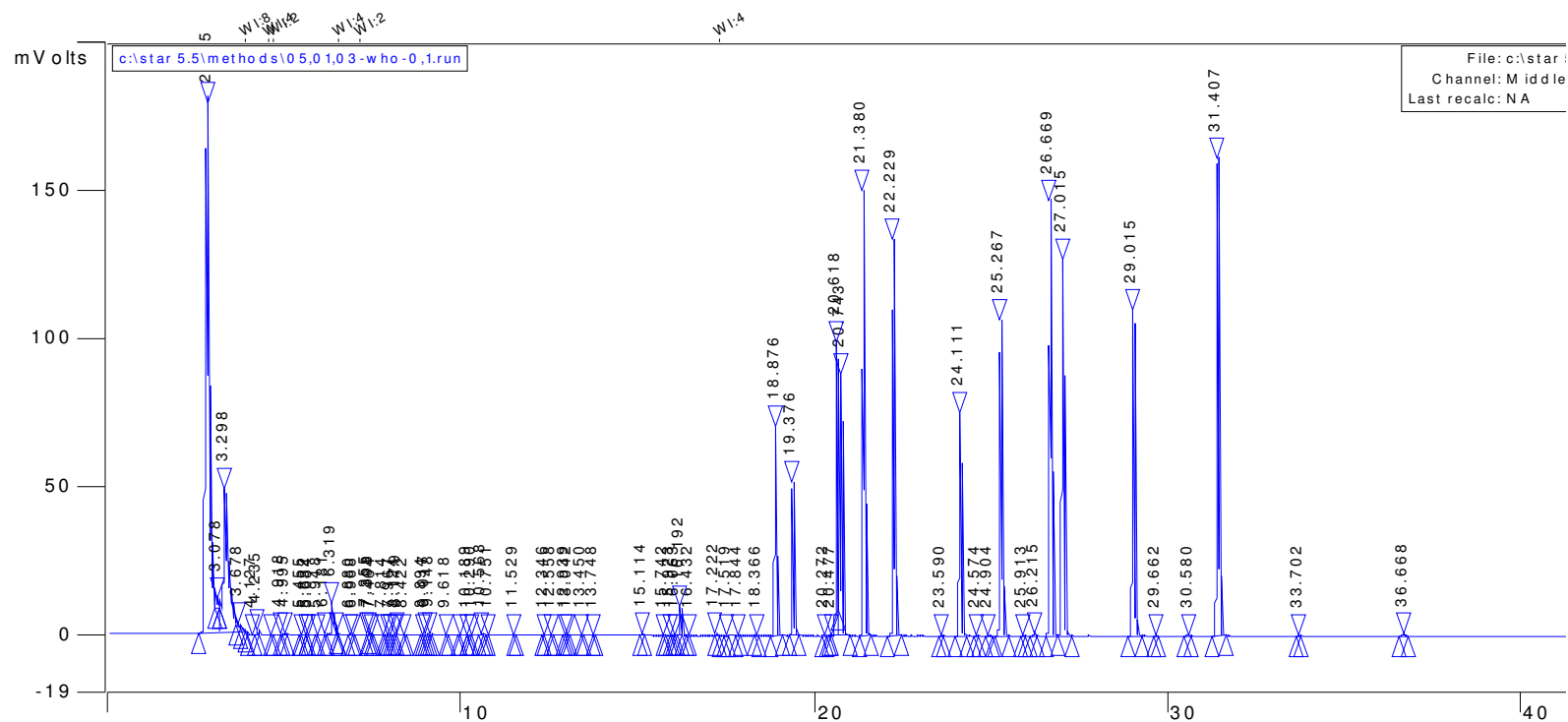


PHB standarts: WHO Congener Mix C-WHO-01 with 12 PCB. Vielu izdalīšanas laiki: PHB81 – 18.88min, PHB 77 – 19.32 min, PHB 123 – 20.62 min, PHB 118 – 20.82 min, PHB 114 – 21.38 min, PHB 105 – 22.23 min, PHB126 – 24.11 min, PHB 167 – 25.3 min, PHB 156 - 26.72 min, PHB 157 – 27.02 min.

9. pielikums

Putekļu ekstrakta hromatogramma ar masspektrometru salīdzinot ar polihlorēto bifenilu standarta hromatogrammu

PHB 12 komponenti

10_03_09_PHB_12_1

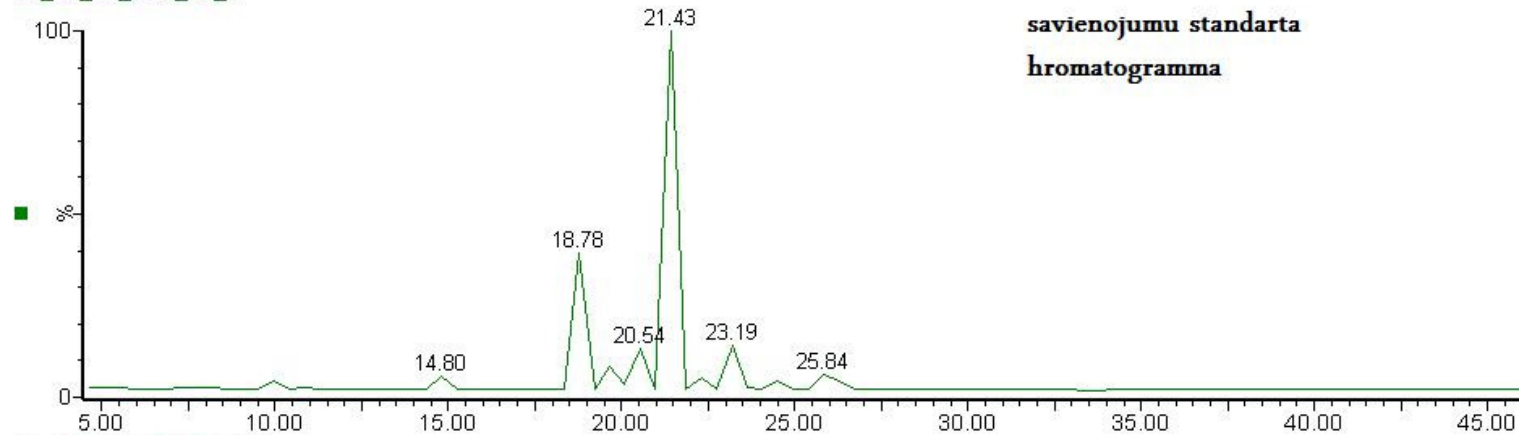

10_03_09_ekstrakts_25

